Total Anchorage Devices (TADs)

ORTHODONTIC MINI-SCREWS

Mir	ni-Sc	FOW	c		
1.111	11-30		2		

38-39

39

Tools



TRANSMUCOSAL SECTION



Scan the QR Code to visit GDT Implants website





Mini Screws

DENTAL IMPLANTS

Mini Screws

ITS



Orthodontic Mini Screws (TADs)

TEMPORARY ANCHORAGE DEVICES

Orthodontic Mini Screws - Temporary Anchorage Devices (TADs) are used to facilitate the movement of teeth through what is called skeletal Anchorage.

TADs Appliances are placed directly into the jawbone then teeth are pulled towards the TADs via elastic or spring traction.

When installed, the Mini screw is in constant and direct contact with the bone and does not possess a periodontal ligament.

As a result, it does not move when orthodontic force is applied, and therefore can be used as "absolute anchorage".

MATERIAL: Ti-6AL-4VELI (Titanium grade 5).

PACKEGING: Individual sterile tube, in a double blister package.


SPECIFICATIONS:



HEXAGONAL HEAD WITH A HOLE

• Hole 0.8mm (suitable for orthodontic wire up to .021 "x.025") for fixation and traction.

TRANSMUCOSAL SECTION

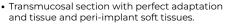

TITANIUM GRADE 5 ALLOY

• Combines properties such as low specific weight, high strength to weight ratio, very high corrosion resistance, and excellent general biocompatibility.

APICAL

- Active point with self-piercing tip profile, specially designed to offer stability secure fit in bone texture and resistance to applications of immediate orthodontic loads.
- Self-tapping.
- Self-drilling.





Mini Screws	DENI			Dental Implants	
Ø <u>2.2mm</u>	ТҮРЕ	TS (mm)	D (mm)	L (mm)	CODE
ل	Chart Namb	1	Ø1.5	6	ORT106
$\overline{\mathbf{O}}$	Short Neck			8	ORT108
	Standard	2	Ø1.5	6	ORT206
Transmucosal Section				8	ORT208
				10	ORT210
E				12	ORT212
H				14	ORT214
F	Wide	2	Ø2	10	ORT410
				12	ORT412
Ft D	Cross-Head	1	Ø1.5	6	ORT316
				8	ORT318
H	Buccal Shelf	4	Ø2	10	ORT215
7 ·				12	ORT216
	Zygomatic	8	Ø2	9	ORT217
HEXAGON	NAL KEY SHO	ORT/LC	DNG	L (mm)	CODE
	25	ORTHKS			
	60	ORTHKL			
SPEAF	L (mm)	CODE			
	60	ORTSTP			
HA	L (mm)	CODE			
	25	ORTHFG			
DRIVER H	L (mm)	CODE			
	93	ORTHTL			
	CODE				

ORTHODONTIC KIT

CODE - ORTKIT



The kit is made of engineering polymer with high resistance to temperatures of high autoclaves temperatures. The tools have standard coupling allowing docking in reducing surgical against angles.

The key to fitting is equipped with contain Nitinol o´ring, which facilitates engagement and disengagement of Miniscrew and prevents the wear of material sterilizations.